



Järven Rescue Sheet "Fast evacuation saves lives"


The rescue sheet, which is always in place between the mattress and the bed, guarantees fast and safe evacuation of bedridden patients in case of fire or other danger. A single person is all that's needed to evacuate a patient. The rescuer simply grabs the handles and can pull the patient to safety.

The Järven Rescue Sheet meets the only standard in the world for evacuation sheets (SS 876 00 39). This standard guarantees that the handles can take a load of 1000 Newton and that the sheet glides easily across all types of surfaces.


The most sold rescue sheet on the market.

More than 100,000 Järven rescue sheets have been sold around the world. Our rescue sheet has been used in a number of successful evacuations where personnel have quickly and effectively saved patients from life-threatening situations.

The Järven Rescue Sheet is produced with a high quality fabric that glides as well on wet as on dry surfaces.

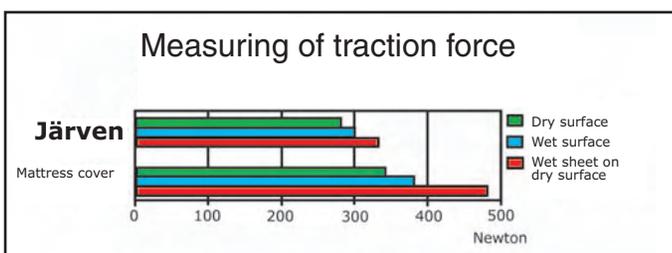

SAFETY - STRENGTH - HYGIENE

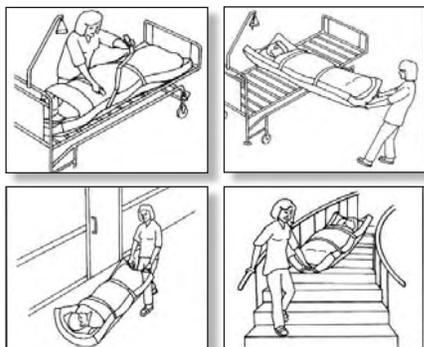

Trade Mark protection No



Always use the Järven Rescue Sheet between the mattress and bed. It's a guarantee for quick evacuations of patients in case of fire or other danger.

The Järven Rescue Sheet has an extra long handle at each short end that allows personnel to evacuate with their back straight, reducing the physical strain.

PRODUCT INFORMATION

- Symmetrical design allows evacuation from both foot and head end.
- Does not interfere with the daily bed making.
- Glides just as well on wet as on dry surfaces.
- X-ray transparent.
- Sturdy pull handles.
- Handles for moving laterally and for daily handling of the mattress.
- Velcro straps for fastening the patient securely during the evacuation.
- Does not discolour hygiene mattresses.
- Flame resistant in accordance with SS 876 00 01.
- Meets the only standard in the world for rescue sheets (SS 876 00 39).

WHEN CLEANING THE RESCUE SHEET IN BED

Use a soap-based cleaning product and dilute in accordance with manufacturer's instructions. Never use concentrate. Cleaners with chlorine may not be used. Disinfect using alcohol-based cleaners with surfactants.

WASHING RECOMMENDATIONS

- Wash separately or in a washing bag at max 90° C.
- Tumble dry at max 70° C.
- Softeners may not be used.
- **NOTE:** The Velcro straps must be secured together before washing the sheet, otherwise the Velcro could damage the sheet and other fabrics and the lint from other fabrics could fasten in the Velcro. This could diminish the Velcro's ability to secure patients during an evacuation.

Art.No.	For mattress width
9900	78-84 cm
9901	70-77 cm
9902	85-90 cm

For other sizes, please indicate measurements when ordering.



(Registered design no. 60487)

Järven's evacuation video and DVD are available for purchase.